

# BOA DUEL TLSO

THORACIC-LUMBO-SACRAL ORTHOSIS BRACE

DONJOY®




# BOA DUEL TLSO

## Features and Benefits

- Rigid Anterior and Posterior panels:
  - Restrict forward flexion and trunk rotation
  - Can be removed and thermoformed for customized fit
- Telescoping Sternal Y-Bar can be:
  - Easily adjusted for length
  - Bent to desired angle for optimal compression
  - Locked into desired fitting for patient safety and compliance
- Easy to use Features:
  - Quick-release shoulder buckles combined with dual action pull handles make for easy donning and doffing
  - Dual compression panels for secure fit and optimal circumferential compression
  - Simple hook-loop closures
- Elastic shoulder straps provide comfort and prevent the brace from “riding up” when sitting
- Lightweight, breathable spacer fabric liner provides comfortable fit
- TLSO comes in 5 sizes, scaled to a single waist size measurement for easy and accurate fitting
- Lifetime guarantee on the proprietary compression panels

Recommended HCPCS Code: L0460

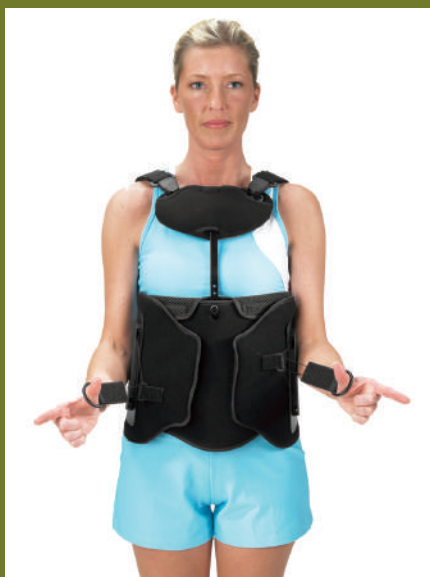

## BOA DUEL TLSO (Thoracic–Lumbo–Sacral Orthosis) Brace

The **Duel TLSO** is a unique prefabricated orthosis (patent pending) for post-operative, post-traumatic, or chronic spinal pathology. The Duel TLSO provides full-circumferential rigid support from the scapular thoracic region of the spine down to the Sacroccygeal junction. The Duel incorporates a telescoping sternal Y-bar that prevents excessive movement. The design works to create a three-point force system consisting of posteriorly directed force from the sternal and pubic pads; and an anteriorly directed force from the thoracolumbar pad. The pressure restricts forward flexion and encourages a hyperextension posture. The posture increases lumbar lordosis and stabilizes or locks the spine to resist lateral and rotary movements.



BOA Duel's Unique Dual Compression Action

## INDICATIONS FOR USE

- Post-Operative Fusion
- Post-Operative Laminectomy
- Post-Operative Disectomy
- Compression Fractures
- Degenerative Disc Disease
- Scoliosis
- Spondylolisthesis
- Spondylolysis
- Spinal Stenosis
- IDET Procedure
- Osteoporosis
- Chronic Low Back Pain

## BOA DUEL TLSO


### Ordering Information

Part Number	Description	Size
11-2291-X	BOA DUEL TLSO	S–XXL

### Sizing Chart

Size	Waist	Size Code (x)
Small	25" – 30"	2
Medium	30" – 35"	3
Large	35" – 42"	4
X-Large	42" – 48"	5
XX-Large	48" – 57"	6

## BOA FAMILY OF PRODUCTS

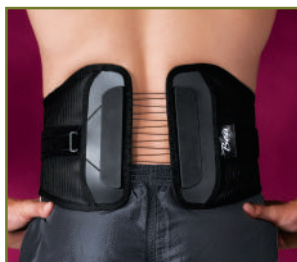

BOA Lite



BOA 8" & 10"



BOA with Chairback



djortho 1430 Decision Street, Vista, CA 92081-8553  
 t800.336.6569 f800.936.6569  
 www.donjoy.com © 2006 dj Orthopedics, LLC



It is within the sole discretion of the customer to determine the appropriate billing code for a product, as well as, whether the use of a product complies with medical necessity and other documentation requirements of the payor. dj Orthopedics, LLC accepts no responsibility whatsoever in this regard. 00-2171 REV A 361